


EPHB1 TE Assay

Assay Format: ADH (1)
NanoBRET™ Tracer: K-9 (2)
100X [Tracer]: 66 μ M in DMSO
Final [Tracer]: 0.66 μ M
Assay Category: High Window (3)
Z': 0.95


Materials Needed

- EPHB1-NanoLuc® Fusion Vector (Cat.# NV3071)
- Transfection Carrier DNA (kit component; 2)
- NanoBRET™ TE Intracellular Kinase Assay, K-9 (Cat.# N2630, N2631, N2830; 2)


Overview of the NanoBRET TE Assay

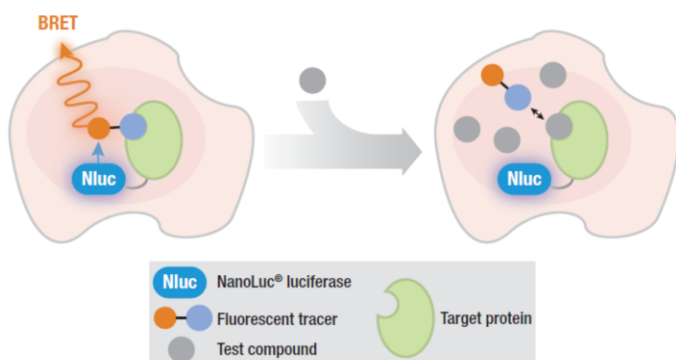

EPHB1-NanoLuc Tracer K-9 Dose Response

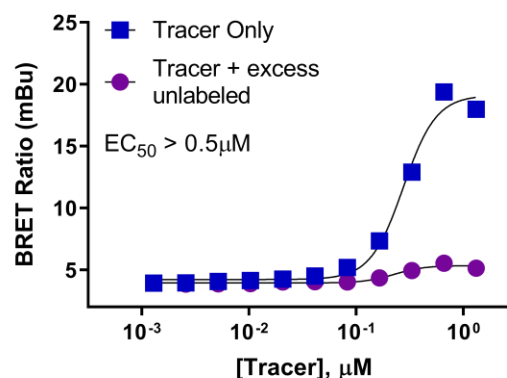


EPHB1-NanoLuc Optimizing [Tracer K-9]

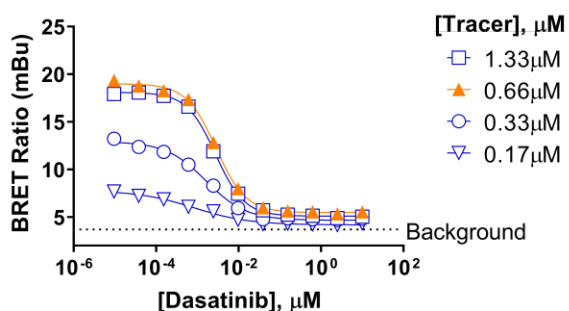

	1.33 μ M	0.66 μ M	0.33 μ M	0.17 μ M
IC ₅₀	0.0028	0.0028	0.0016	0.00083

EPHB1-NanoLuc Example Compound Profiling

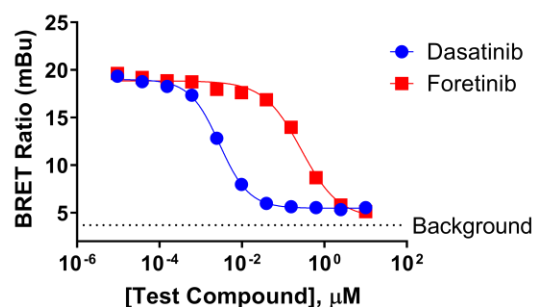

	Dasatinib	Foretinib
IC ₅₀	0.0028	0.27

Example NanoBRET™ Tracer K-9 data in HEK293 cells transiently expressing EPHB1-NanoLuc®

Top Left Panel: Overview of the NanoBRET™ TE Assay. **Other Panels:** HEK293 cells were first transfected with EPHB1-NanoLuc® Fusion Vector and Transfection Carrier DNA (1 to 9 transfection ratio) and then were subsequently resuspended in OptiMEM prior to seeding into 96-well plates: **Top Right Panel:** Tracer affinity was measured by treating the cells with increasing concentrations of tracer in the presence or absence of a molar excess of unlabeled compound. **Bottom Left Panel:** The apparent cellular affinity of the unlabeled compound was measured at multiple fixed concentrations of the tracer, where the IC₅₀ at the recommended tracer concentration is depicted in orange (4). **Bottom Right Panel:** An example compound profiling experiment at the recommended tracer concentration is provided.



NanoBRET™ TE Intracellular Kinase Assay

Notes:

- (1) This assay is run in the adherent format. See the **NanoBRET™ TE Intracellular Kinase Assay, Adherent Format** Technical Manual (#TM598) for protocol details.
- (2) NanoBRET™ Tracer K-9 is supplied within the **NanoBRET™ TE Intracellular Kinase Assay, K-9** products (N2630, N2631, N2830). Additional assay components are supplied within these kits, including the NanoBRET™ Nano-Glo® Substrate, Extracellular NanoLuc® Inhibitor, and TE Tracer Dilution Buffer. Additionally, NanoLuc®-PLK2 Fusion Vector (as a control) and Transfection Carrier DNA are provided in products N2630 and N2631. Transfection carrier DNA products (E4881 and E4882) are also available separately. For full details, please see the Promega website or technical manual for these products.
- (3) Assay category is defined by the assay window at the recommended tracer concentration. It is detailed in section 7.A within the **NanoBRET™ TE Intracellular Kinase Assay, Adherent Format** Technical Manual (#TM598).
- (4) See section 5 of **NanoBRET™ TE Intracellular Kinase Assay, Adherent Format** Technical Manual (#TM598) regarding approaches to achieve quantitative analysis of test compound affinity.
- (5) NanoBRET™ TE Intracellular Assays have also been applied to Residence Time analysis. For a kinase example, please refer to Forster, M. *et al.* For an HDAC example please refer to Robers, M.B. *et al.*

References:

Vasta, J. D. *et al.* (2018) Quantitative, wide-spectrum kinase profiling in live cells for assessing the effect of cellular ATP on target engagement. *Cell Chem. Biol.* **25**, 206.

Robers, M.B. *et al.* (2015) Target engagement and drug residence time can be observed in living cells with BRET. *Nature Comm.* **6**, 10091.

Forster, M. *et al.* (2016) Selective JAK3 inhibitors with a covalent reversible binding mode targeting a new induced fit binding pocket. *Cell Chem. Biol.* **23**, 1335.

This protocol was developed by Promega Scientists and is intended for research use only.

Users are responsible for determining suitability of the protocol for their application.

Further information can be found in Technical Manual #TM598, available at: www.promega.com/protocols