


## MAP4K2 TE Assay

**Assay Format:** ADH (1)  
**NanoBRET™ Tracer:** K-10 (2)  
**100X [Tracer]:** 6.3 μM in DMSO  
**Final [Tracer]:** 0.063 μM  
**Assay Category:** High Window (3)  
**Z':** 0.85


### Materials Needed

- NanoLuc®-MAP4K2 Fusion Vector (Cat.# NV1621)
- Transfection Carrier DNA (kit component; 2)
- NanoBRET™ TE Intracellular Kinase Assay, K-10 (Cat.# N2640, N2641, N2840; 2)


### Overview of the NanoBRET TE Assay

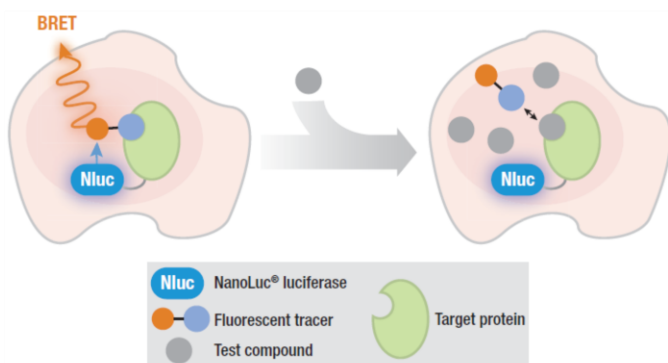

### NanoLuc-MAP4K2 Tracer K-10 Dose Response

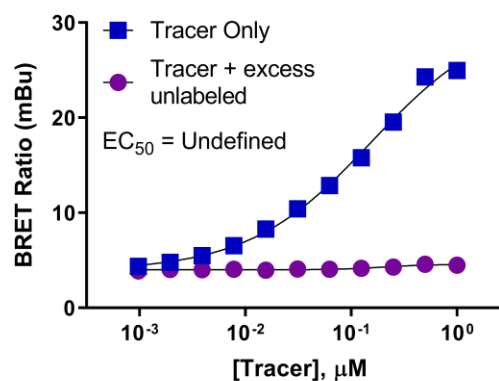


### NanoLuc-MAP4K2 Optimizing [Tracer K-10]

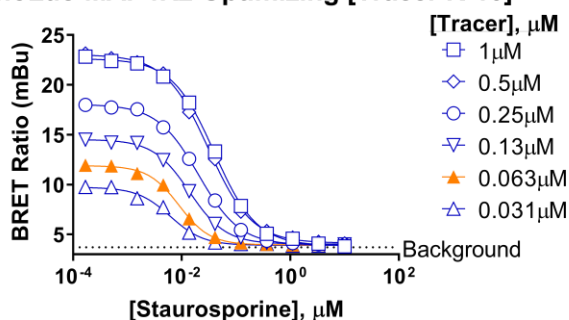

	1 μM	0.5 μM	0.25 μM	0.13 μM	0.063 μM	0.031 μM
IC50	0.040	0.031	0.020	0.014	0.0089	0.0062

### NanoLuc-MAP4K2 Example Compound Profiling

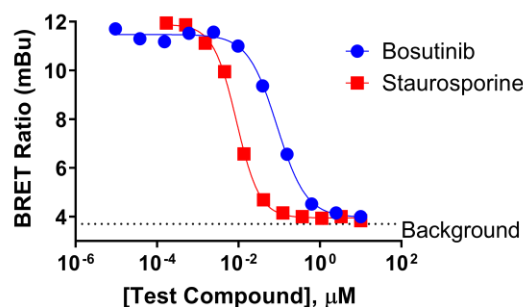

	Bosutinib	Staurosporine
IC50	0.088	0.0089

### Example NanoBRET™ Tracer K-10 data in HEK293 cells transiently expressing NanoLuc®-MAP4K2

**Top Left Panel:** Overview of the NanoBRET™ TE Assay. **Other Panels:** HEK293 cells were first transfected with NanoLuc®-MAP4K2 Fusion Vector and Transfection Carrier DNA (1 to 9 transfection ratio) and then were subsequently resuspended in OptiMEM prior to seeding into 96-well plates: **Top Right Panel:** Tracer affinity was measured by treating the cells with increasing concentrations of tracer in the presence or absence of a molar excess of unlabeled compound. **Bottom Left Panel:** The apparent cellular affinity of the unlabeled compound was measured at multiple fixed concentrations of the tracer, where the IC<sub>50</sub> at the recommended tracer concentration is depicted in orange (4). **Bottom Right Panel:** An example compound profiling experiment at the recommended tracer concentration is provided.



## NanoBRET™ TE Intracellular Kinase Assay

### Notes:

- (1) This assay is run in the adherent format. See the **NanoBRET™ TE Intracellular Kinase Assay, Adherent Format** Technical Manual (#TM598) for protocol details.
- (2) NanoBRET™ Tracer K-10 is supplied within the **NanoBRET™ TE Intracellular Kinase Assay, K-10** products (N2640, N2641, N2840). Additional assay components are supplied within these kits, including the NanoBRET™ Nano-Glo® Substrate, Extracellular NanoLuc® Inhibitor, and TE Tracer Dilution Buffer. Additionally, LIMK1-NanoLuc® Fusion Vector (as a control) and Transfection Carrier DNA are provided in products N2640 and N2641. Transfection carrier DNA products (E4881 and E4882) are also available separately. For full details, please see the Promega website or technical manual for these products.
- (3) Assay category is defined by the assay window at the recommended tracer concentration. It is detailed in section 7.A within the **NanoBRET™ TE Intracellular Kinase Assay, Adherent Format** Technical Manual (#TM598).
- (4) See section 5 of **NanoBRET™ TE Intracellular Kinase Assay, Adherent Format** Technical Manual (#TM598) regarding approaches to achieve quantitative analysis of test compound affinity.
- (5) NanoBRET™ TE Intracellular Assays have also been applied to Residence Time analysis. For a kinase example, please refer to Forster, M. *et al.* For an HDAC example please refer to Robers, M.B. *et al.*

### References:

Vasta, J. D. *et al.* (2018) Quantitative, wide-spectrum kinase profiling in live cells for assessing the effect of cellular ATP on target engagement. *Cell Chem. Biol.* **25**, 206.

Robers, M.B. *et al.* (2015) Target engagement and drug residence time can be observed in living cells with BRET. *Nature Comm.* **6**, 10091.

Forster, M. *et al.* (2016) Selective JAK3 inhibitors with a covalent reversible binding mode targeting a new induced fit binding pocket. *Cell Chem. Biol.* **23**, 1335.

This protocol was developed by Promega Scientists and is intended for research use only.

Users are responsible for determining suitability of the protocol for their application.

Further information can be found in Technical Manual #TM598, available at: [www.promega.com/protocols](http://www.promega.com/protocols)